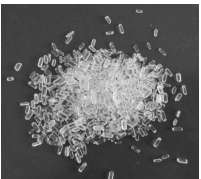




**% Water in Magnesium Sulfate**



$\text{MgSO}_4 \cdot 7 \text{H}_2\text{O}$





$\text{MgSO}_4 + 7 \text{H}_2\text{O}$

Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---

---

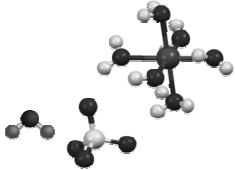
---

---

---

**Hydrates**

**Absorb water from atmosphere**  
**Water becomes associated with structure**



$\text{MgSO}_4 \cdot 7 \text{H}_2\text{O}$

Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---


---

---

---


---

**Hydrates**



**Compound + water → hydrate**  
**reversible**

**Hydrate → compound + water**



$\text{CuSO}_4 \cdot 5 \text{H}_2\text{O} \rightleftharpoons \text{CuSO}_4 + 5 \text{H}_2\text{O}$

$\text{MgSO}_4 \cdot 7 \text{H}_2\text{O} \rightleftharpoons \text{MgSO}_4 + 7 \text{H}_2\text{O}$

**Reactants and Products are chemically different:**  
**Color change indicates chemical change**  
**These reactions represent chemical changes**

Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---

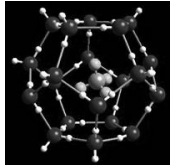
---

---

---

---

## Methane Hydrates (Clathrates)



At cold temperatures:  
Methane trapped by ice  
Abundant in tundra and ocean bottoms  
Global warming releasing the methane



Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

---

---

---

---

## Desiccants

Compounds that absorb water to form hydrates



Used to protect variety of commercial products

Keep desiccants in containers until contents consumed

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

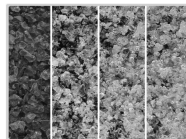
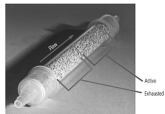
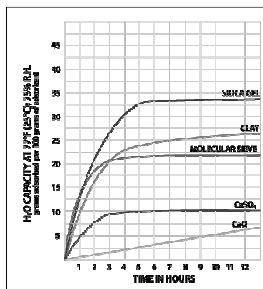
---

---

---

---

## Desiccants Five Common Types



Color Change shows absorption

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

---



---


---

---


## Lab Desiccators


Used to Keep Sensitive Chemicals Dry



The Desiccator





**Often stored under vacuum and sometimes in the cold (-78 °C)**

Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---

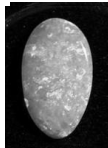

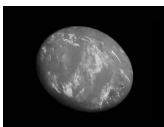
---

---


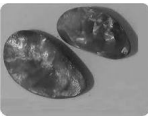
---

---

## Color of Many Gems From Hydrates

**Never store in dehydrating conditions  
Or  
Underwater**

Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---

---


---

---

---

## Composition Calculations


Find the percent water present in the hydrate  $\text{CuSO}_4 \cdot 5 \text{H}_2\text{O}$

Cu = 63.55		
S = 32.07	10 H = 10.08	
4 O = 64.00	5 O = 80.00	
<b>Total = 159.62</b>	<b>Total = 90.08</b>	

$$\% \text{ Water} = \frac{90.08}{159.62 + 90.08} \times 100 = 36.08$$

(249.70)

**36.08 % = 0.3608**



Copyright Larry P. Taylor, Ph.D. All Rights Reserved LPT

---

---

---

---

---

---

---

---

### Composition Calculations

Find the percent water present in the hydrate  $\text{MgSO}_4 \cdot 7 \text{H}_2\text{O}$

$$\text{Mg} = 24.31$$

$$\text{S} = 32.07$$

$$4 \text{ O} = 64.00$$

$$\text{Total} = 120.38$$

$$14 \text{ H} = 14.11$$

$$7 \text{ O} = 112.00$$

$$\text{Total} = 126.11$$



$$\% \text{ Water} = \frac{126.11}{120.38 + 126.11} \times 100 = 51.24$$

(246.49)



$$51.24 \% = 0.5124$$

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

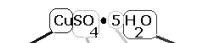
---

---


---

### Naming Hydrates


Anhydrous (without water) name “• n H<sub>2</sub>O’s”




copper (II) sulfate penta hydrate

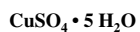

copper (II) sulfate pentahydrate

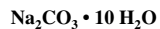

sodium carbonate decahydrate

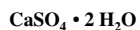

calcium sulfate dihydrate

- Indicates distinct chemical entities held together



Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

Desiccants keep my wraps dry



Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---




---

---

---

---

## % Water in Magnesium Sulfate Heptahydrate

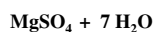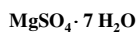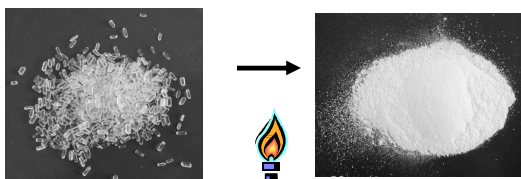

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

---

---

## Hydrate Lab

### Purpose

Determine the percentage of water in a given hydrate

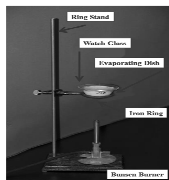
### Procedure



Weigh materials "by difference"

Weight of evaporating dish, watch glass, & hydrate  
- Weight of evaporating dish & watch glass

Weight of hydrate



Water driven away by heat  
Watch glass minimizes splattering  
Heat until all water is gone



Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

---

---

## Calculations

Show all your work

Mass of a any substance (weighing by difference)

(Substance g + Container g) - Container g = Substance g

Theoretical Water loss: initial heptahydrate x % H<sub>2</sub>O

Water Lost: Initial - final weight of the magnesium hydrate

Experimental % Water:  $\frac{(\text{mass H}_2\text{O lost})}{(\text{mass initial heptahydrate})} \times 100$

## Results

Tabulate the answers to your calculations

## Conclusion

State % water in MgSO<sub>4</sub> · 7 H<sub>2</sub>O

Compare your experimental value to the theoretical

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

---

---

## Determine n



Calculate the value of n for  $\text{MgSO}_4 \cdot n \text{H}_2\text{O}$

N is the ratio of moles water to moles anhydrous salt  
Experiment measures grams ... need moles for this ratio

Convert grams water lost to moles (via molar mass of one  $\text{H}_2\text{O}$ )  
Convert grams  $\text{MgSO}_4$  remaining to moles (via molar mass  $\text{MgSO}_4$ )

$$n = \frac{\text{Moles water lost}}{\text{Moles anhydrous magnesium heptahydrate}}$$

n is closest small, whole number

Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---

## Let's Boldly Go Explore Today's Lab



Copyright Larry P. Taylor, Ph.D. All Rights Reserved

LPT

---

---

---

---

---

---

---

---